



Oclean

Product User Manual

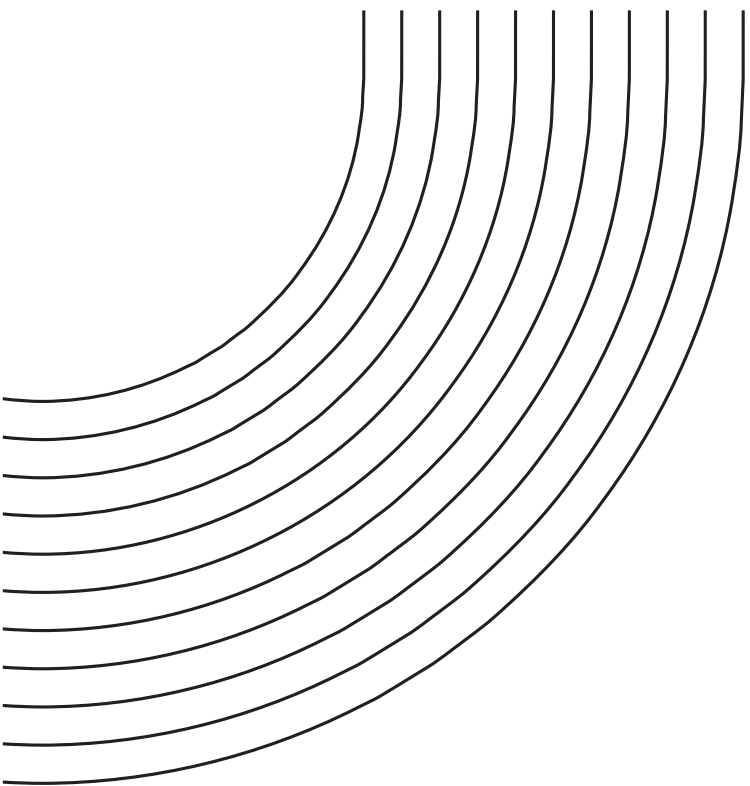


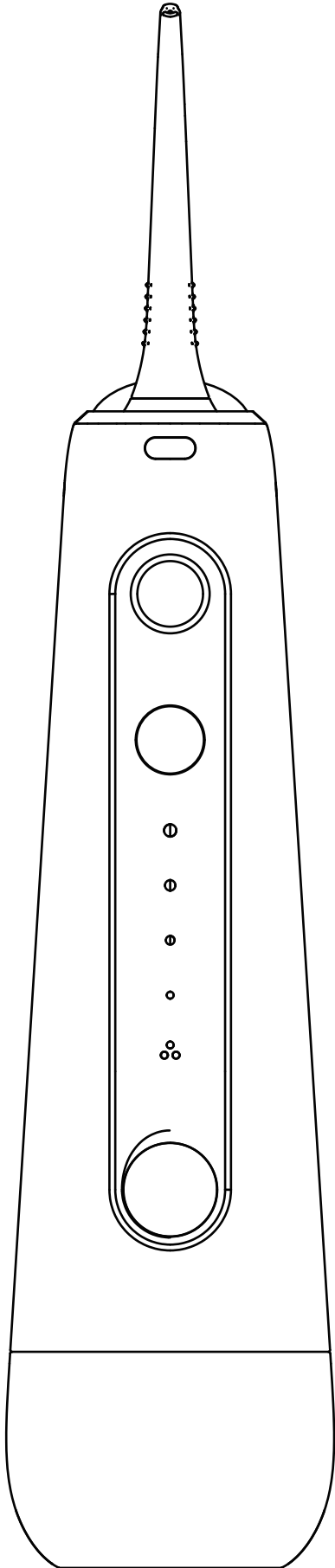

Water Flosser



Standard Nozzle

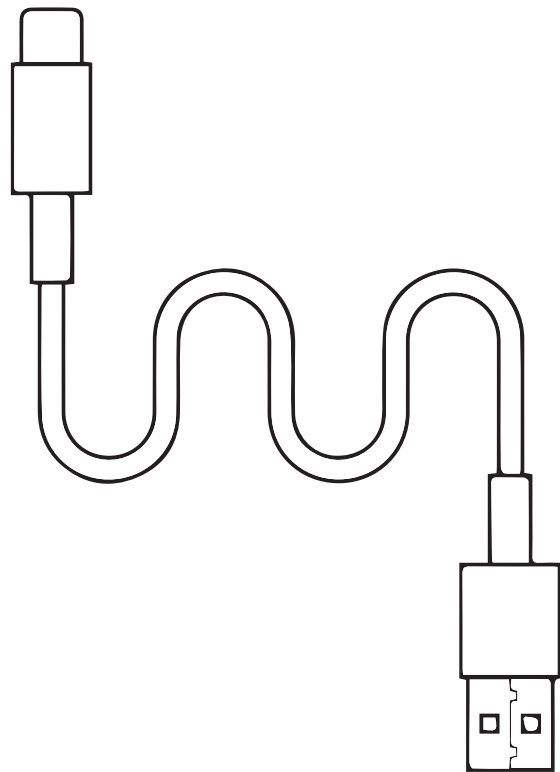

Charging Cable

Thank you for using Oclean oral irrigator.

Oclean W10 Water Flosser

Important Safety Precautions

1. This product can only be used with Oclean-specific nozzles. Oclean is not responsible for any damage or loss caused by the use of any third-party accessories.
2. Do not use this product without a nozzle installed.
3. It is not recommended for children under 12 years old. Children should not play with this device.
4. Patients who have just undergone oral surgery should consult a dentist before using this product.
5. If the gums continue to bleed after using this product for more than 2 weeks, please consult a dentist.
6. Do not place the tip of the nozzle close to the gums or teeth to avoid injury.
7. Do not use it for purposes other than oral cleaning, otherwise it may cause accidental injury when water sprays into the eyes, nose,

ears or throat.

8. Do not use iodine, saline or concentrated essential oils that are not soluble in water in your product. Using these liquids will reduce the performance and shorten the service life of this product. Do not inject highly concentrated mouthwash into this product. If you need to use mouthwash, please dilute it in proportion before feeding into this product. After using the diluted mouthwash, you need to rinse it again with clean water to clean the inner tube.

9. Do not immerse the oral irrigator in water. Please wipe the device with a dry cloth after use.

10. If the product is not to be used for a long time, the water tank should be emptied and dried, the nozzle, body and charging port should be cleaned, and it should be stored in a cool and dry place.

11. Do not disassemble this product or replace its parts to avoid damage or danger.

12. The new device may have a slight plastic smell in the first few uses, which is normal and the smell will disappear after a period of use. If the product you purchased is faulty or defective, please contact our service hotline: service@o-clean.com

Flash in red:
Low battery alarm

Flash in white:
Charging status

Keep on in white:
Battery is fully charged

Intensive Mode:
thorough and deep
clean intensity

Standard Mode:
for daily general clean
intensity

Gentle Mode:
for those of sensitive
gums brushing intensit

Pulsating Mode:
rhythmed gum
massage intensity


On-demand Mode:
floss just the way
you want intensity

**Nozzle installation
button:**
Press the button and
insert the nozzle into
the oral irrigator body.
Press the button again
to remove the nozzle
from the oral irrigator
body.


Power button:
Press the Power
button to start the oral
irrigator, and the white
indicator will stay on.
Press the Power
button again to turn
off the oral irrigator.

Mode button


**Type-C Charging
Port:**
Open the cover and
connect the USB
charging cable to
start charging.

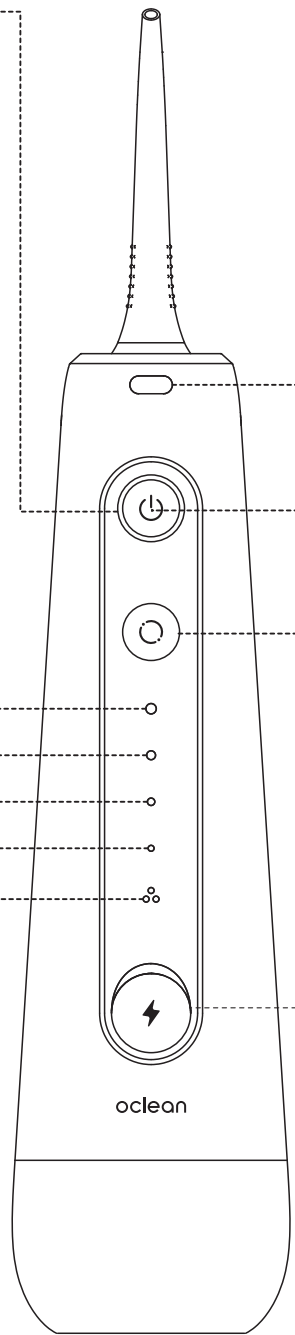


Instruction For Use

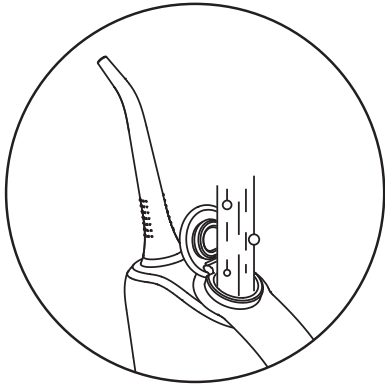

Open the cover and fill the water tank with clean water (recommended water temperature below 40°C), then close the cover.

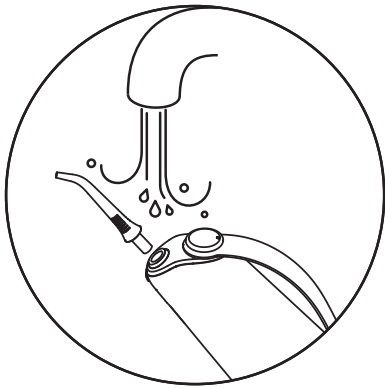

Remove the nozzle from the device, clean it, and remove any residue after each use.

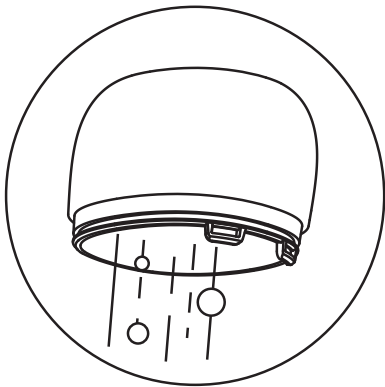

After use, press the button on the bottom cover of the water tank and slowly pull it down to empty the remaining water.

Water filling operation: Press the "Power" button to start using this product. This product has an automatic shutdown function. Ensure there is still water in the tank and press the "Power" button again. To pause during use, press the "Power" button. Press the "Power" button again to continue.

Product Specifications

Product Name	Oclean W10 Portable Oral Irrigator
Model	W6010
Battery Capacity	1400mAh(**)Lithium ion battery
Battery Life	30Days(*)
Water Compartment Capacity	200ml(±5%)
Charging Time	5Hours(***)
Product Size	72×55×267mm(nozzle height included)
Input	DC 5V=1A
Charging	USB Type-C charging cable
Product Net Weight	278g

*Note:Use the irrigator once a day in impulse mode, using one compartment of water for each use.

Theoretically, the battery capacity has a margin of ±50mAh.

Use a standard adapter of 5V/1A.

Battery Removal Instructions

1. Ensure the product is powered off before removing the battery.
2. Dispose of the battery safely after removal.
3. Battery removal should be performed by a professional following these steps:
 - Remove the water tank, nozzle, and hand strap from the bottom of the machine.
 - Use a screwdriver to pry open the top decorative piece from the water inlet, then unscrew the fastening screws from the top with the screwdriver and remove the top fixture of the machine.
 - Use pliers to cut open the latch at the bottom of the shell, then use a screwdriver to push out the core assembly from the top water outlet.
 - Use pliers to cut the positive and negative battery leads, then use a screwdriver to loosen the battery holder screws and remove the battery.